

Stereoscope Instructions

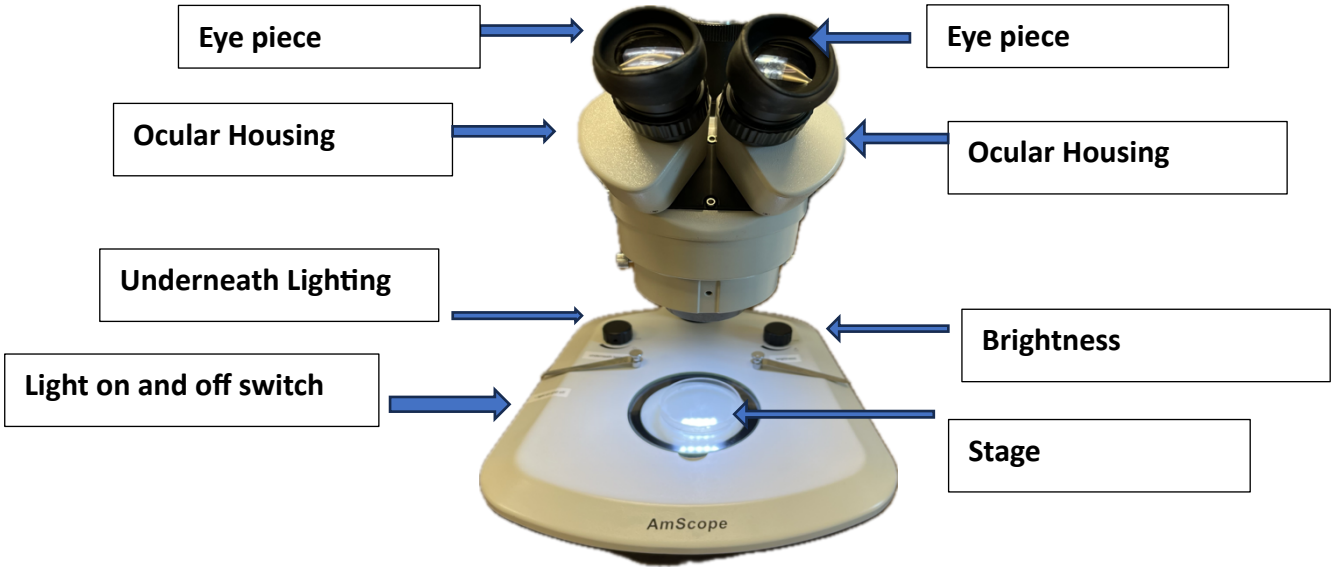

- ☐ Place Specimen on a plastic tray in the center of the **stage,** not directly on the glass.
- ☐ Turn on using the **light on and off** switch, this turns on top lighting. Top lighting is ideal for solid specimens whilst **underneath lighting** maybe preferable for translucent samples. The **brightness** can also be adjusted.
- ☐ Before looking through the eye pieces make sure the lens is at its lowest position using **ONE focusing knob** only**.** **<u>Warning - using more than one focusing knob at once may damage the scope.</u>**
- ☐ With both eyes open look into the eye pieces. Adjust the eye distance by slowly rotating the ocular housings in either direction until one circle of light can be seen.
- ☐ Close one eye. Using the focusing knob adjust until the specimen is in focus, close the other eye and use eye piece to focus.
- ☐ Open both eyes and see your specimen. To zoom in use the zoom knob.
- ☐ After use switch off lights and re-cover the stereoscope.

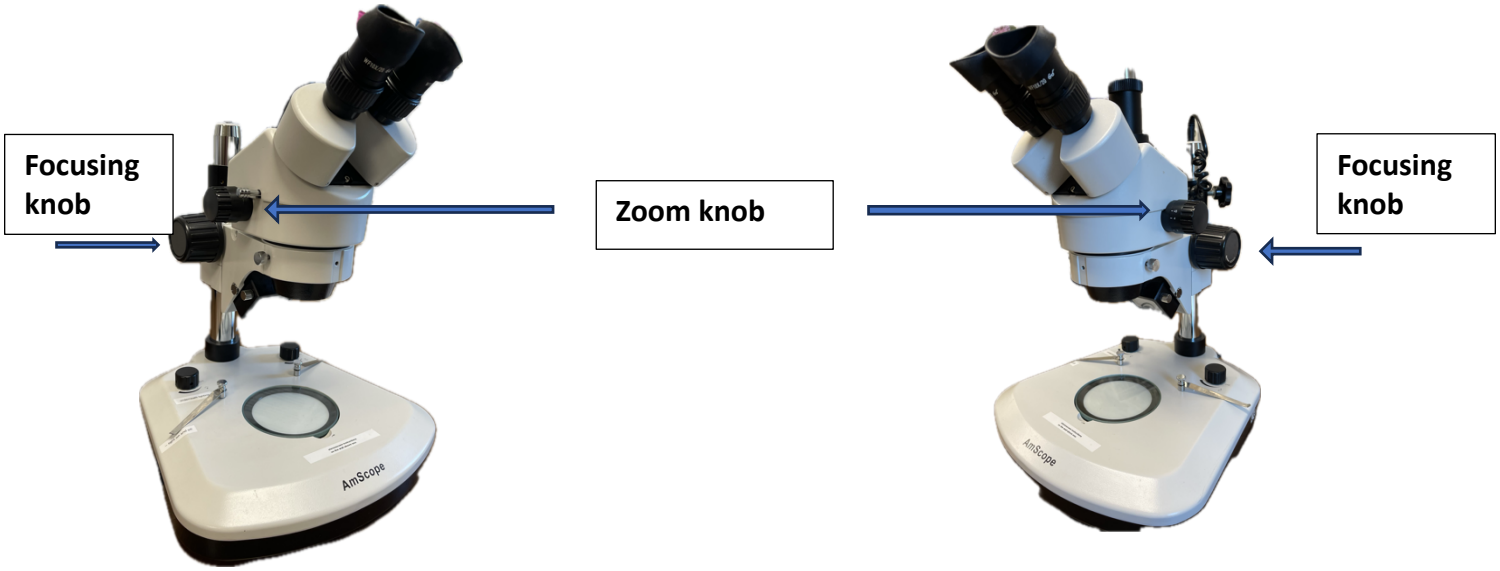

HELPLINE HINT 4